

Product Sheet

H_GLP1R GIPR GCGR Triple Reporter 293 Cell Line


Catalog number: GM-C44406

Version 3.3.1.260804

GIPR, GCGR, and GLP-1R are all Class B GPCRs. Their activation primarily triggers the cAMP-PKA signaling pathway, but they differ in distribution and function, GIPR is mainly in pancreatic β cells and regulates insulin-related functions. GCGR is expressed in the liver and kidneys, promoting glycogen breakdown. GLP-1R is found in pancreatic cells and brain neurons, enhancing insulin release and appetite control. It is a key target for diabetes drugs.

Although traditional single-target therapies have demonstrated certain therapeutic effects, they still struggle to systematically correct the overall dysregulation of metabolic networks. In recent years, dual-target drugs targeting GLP1R and GIPR have developed rapidly. Owing to their significant advantages in "synergistic efficacy and reduced side effects," they have become a key trend in drug development in this field..

The H_GLP1R GIPR GCGR Triple Reporter 293 Cell Line is a clonal, stable cell line that constitutively expresses human GLP1R, human GIPR and human GCGR, along with a signal-dependent expression of a luciferase reporter gene. The binding of Retatrutide to the three receptors activate downstream reporter genes, leading to luciferase expression. The luciferase readout indicates the activation level of the signaling pathway and can thus be used to evaluate the in vitro effects of drugs related to the GLP1R, GIPR and GCGR receptors.

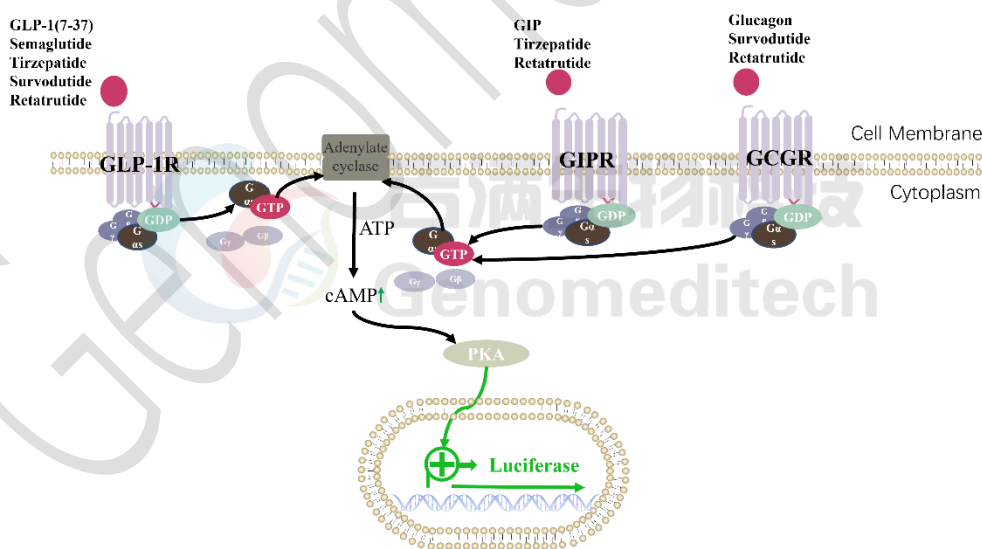

Specifications

Quantity	5E6 Cells per vial, 1 mL
Product Format	1 vial of frozen cells
Shipping	Shipped on dry ice
Storage Conditions	Liquid nitrogen immediately upon receipt
Recovery Medium	DMEM+10% FBS+1% P.S
Growth medium	DMEM+10% FBS+1% P.S+4 µg/mL Blasticidin+400 µg/mL G418+125 µg/mL Hygromycin+0.75 µg/mL Puromycin
Note	None
Freezing Medium	90% FBS+10% DMSO
Growth properties	Adherent
Growth Conditions	37°C, 5% CO ₂
Mycoplasma Testing	The cell line has been screened to confirm the absence of Mycoplasma species.
Safety considerations	Biosafety Level 2
Note	It is recommended to expand the cell culture and store a minimum of 10 vials at an early passage for potential future use.

Materials

Reagent	Manufacturer/Catalogue No.
Puromycin	Genomeditech/ GM-040401
Blasticidin	Genomeditech/ GM-040404
G418	Genomeditech/ GM-040402
Hygromycin	Genomeditech/ GM-040403
Pen/Strep	Thermo/15140-122
Fetal Bovine Serum	ExCell/FSP500
DMEM	Gibco/C11995500BT
Anti-GLP1R hIgG1 Antibody(mAb-36986)	Genomeditech/ GM-51168AB
Anti-H_GIPR hIgG1 Antibody(AMG-133)	Genomeditech/ GM-84915AB
Anti-H_GCGR hIgG2 Antibody(volagidemab)	Genomeditech/ GM-84555AB
GMOne-Step 2.0 Luciferase Reporter Gene Assay Kit	Genomeditech/ GM-040513
Retatrutide	Targetmol/T76279
Glucagon (1-29), bovine, human	MCE/HY-P0082
Tirzepatide	Aladdin/T412109

GIP(Human)

Genscript/RP10795

Semaglutide

aladdin/S774153

Survodutide

MCE/HY-P4146

Genomeditech

Figures

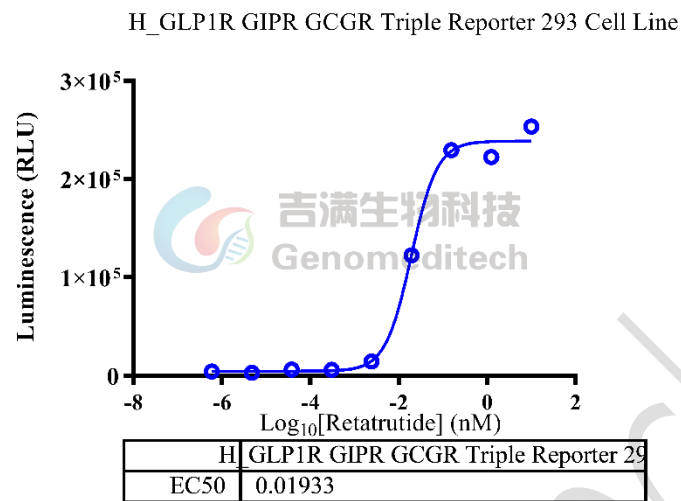


Figure 1 | Response to Retatrutide. The H_GLP1R GIPR GCGR Triple Reporter 293 Cell Line (Cat. GM-C44406) at a concentration of 1.5E4 cells/well (96-well format) was stimulated with serial dilutions of Retatrutide (targetmol/T76279) in assay buffer (DMEM + 1% FBS + 1% P.S) for 6 hours. The firefly luciferase activity was measured using the Luciferase Reporter Assay Kit (Genomeditech). The maximum induction fold was approximately [76.9]. Data are shown by drug molar concentration.

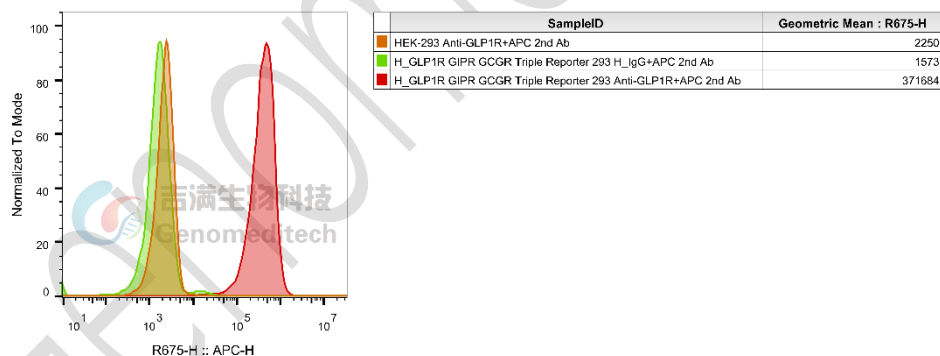


Figure 2 | H_GLP1R GIPR GCGR Triple Reporter 293 Cell Line (Cat. GM-C44406) was determined by flow cytometry using Anti-GLP1R hIgG1 Antibody(mAb-36986) (Cat. [GM-51168AB](#)).

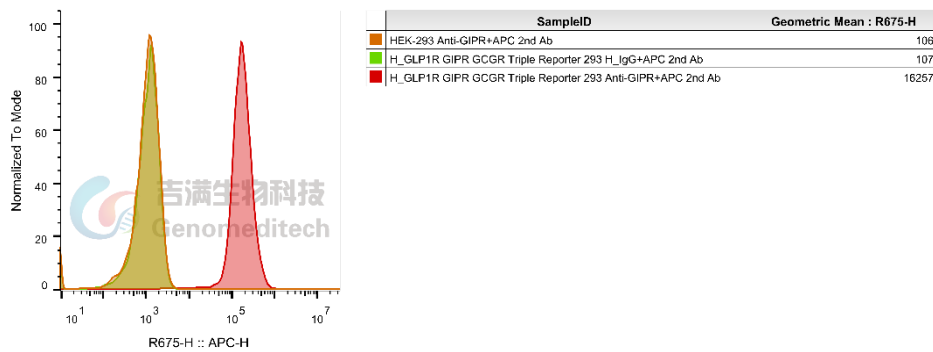


Figure 3 | H_GLP1R GIPR GCGR Triple Reporter 293 Cell Line (Cat. GM-C44406) was determined by flow cytometry using Anti-H_GIPR hIgG1 Antibody(AMG-133) (Cat. [GM-84915AB](#)).

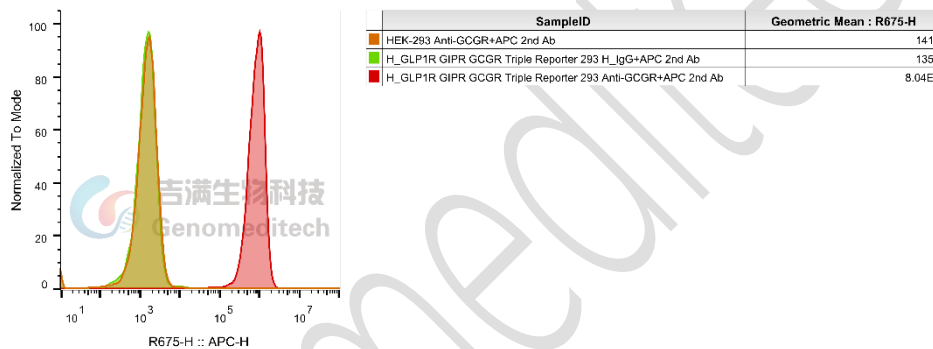


Figure 4 | H_GLP1R GIPR GCGR Triple Reporter 293 Cell Line (Cat. GM-C44406) was determined by flow cytometry using Anti-H_GCGR hIgG2 Antibody(volagidemab) (Cat. [GM-84555AB](#)).

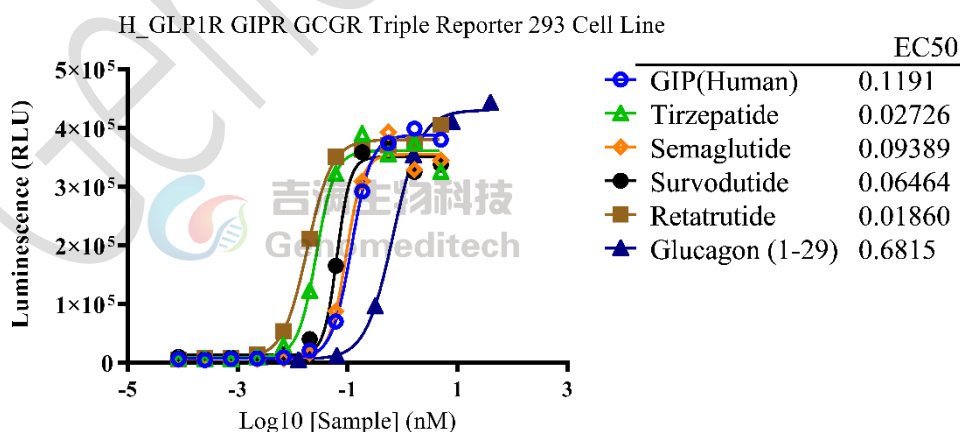


Figure 5 | Response to Glucagon (1-29),GIP, Tirzepatide, Semaglutide, Survodutide, Retatrutide.The H_GLP1R GIPR GCGR Triple Reporter 293 Cell Line (Cat. GM-C44406) at a concentration of 1.5E4 cells/well (96-well format) were stimulated with serial dilutions of Glucagon (1-29),GIP, Tirzepatide, Semaglutide, Survodutide, Retatrutide in assay buffer (DMEM + 1% FBS + 1% P.S) for 6 hours. The firefly luciferase activity was measured using the Luciferase

Reporter Assay Kit (Genomeditech). The maximum induction fold was approximately [63.3],[56.5],[48.8],[51.3],[50.3] and [60.3]. Data are shown by drug molar concentration.

Cell Recovery

Recovery Medium: DMEM+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 - 3 minutes).
- Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- Centrifuge at 176 x g for 3 minutes to collect cells.
- Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- Aliquot 1 mL into each vial.
- Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: DMEM+10% FBS+1% P.S+4 µg/mL Blasticidin+400 µg/mL G418+125 µg/mL Hygromycin+0.75 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- Subculturing is necessary when the cell density reaches 80%. It is recommended to perform subculturing at a ratio of 1:3 to 1:4 every 2-3 days. Ensure that the density does not exceed 80%, as overcrowding can lead to reduced viability due to compression.
- Remove and discard culture medium.
- Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.

- d) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 30 to 60 seconds at 37°C).
- e) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- f) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- g) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- h) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:4 is recommended

Medium Renewal: Every 2 to 3 days

Notes

- a) Upon initial thawing, a higher number of dead cells is observed, which is a normal phenomenon. Significant improvement is seen after adaptation. Once the cells reach a stable state, the number of dead cells decreases after subculturing and the cell growth rate becomes stable.
- b) Ensure that the cell density does not exceed 80%, as overcrowding may lead to reduced viability due to compression.

Related Products

GCGR	
H_GCGR Reporter CHO-K1 Cell Line	H_GCGR Reporter HEK-293 Cell Line
H_GCGR Reporter HEK-293 DDX35TM Cell Line	Cynomolgus_GCGR HEK-293 Cell Line
H_GCGR CHO-K1 Cell Line	H_GCGR HEK-293 Cell Line
Mouse_GCGR HEK-293 Cell Line	
Anti-H_GCGR hIgG2 Antibody(volagidemab)	
Human GCGR(26-136) Protein; His Tag	
GLP1R	
H_GLP1R GIPR Reporter CHO-K1 Cell Line	H_GLP1R Reporter CHO-K1 Cell Line
H_GLP1R Reporter HEK-293 Cell Line	H_GLP1R Reporter HEK-293 DDX35TM Cell Line
H_GLP1R β-Arrestin Reporter CHO-K1 Cell Line	Cynomolgus_GLP1R GIPR CHO-K1 Cell Line
Cynomolgus_GLP1R HEK-293 Cell Line	H_GLP1R CHO-K1 Cell Line
H_GLP1R GIPR CHO-K1 Cell Line	H_GLP1R HEK-293 Cell Line
Mouse_GLP1R GIPR CHO-K1 Cell Line	Mouse_GLP1R HEK-293 Cell Line
Rat_GLP1R HEK-293 Cell Line	
Anti-GLP1R hIgG1 Antibody(mAb-36986)	Anti-H_GLP1R hIgG1 Antibody(glutazumab)
FGFR1	
H_FGF21 Reporter HEK-293 Cell Line	Mouse_FGF21 Reporter HEK-293 Cell Line
Human FGF-21 Protein; His Tag	
CALCA(CGRP):CALCRL RAMP	
H_CALCRL RAMP1 Reporter HEK-293 Cell Line	H_CALCRL RAMP1 Reporter HEK-293 DDX35TM Cell Line

Cynomolgus_CALCRL RAMP1 HEK-293 Cell Line	H_CALCRL RAMP1 CHO-K1 Cell Line
H_CALCRL RAMP1 HEK-293 Cell Line	
Anti-CALCRL RAMP1 hIgG2 Antibody(Erenumab)	
GPR75	
H_GPR75 HEK-293 Cell Line	
GIPR	
H_GIPR Reporter CHO-K1 Cell Line	H_GIPR Reporter HEK-293 Cell Line
H_GIPR Reporter HEK-293 DDX35TM Cell Line	Cynomolgus_GIPR CHO-K1 Cell Line
Cynomolgus_GIPR HEK-293 Cell Line	H_GIPR CHO-K1 Cell Line
H_GIPR HEK-293 Cell Line	Mouse_GIPR CHO-K1 Cell Line
Mouse_GIPR HEK-293 Cell Line	
Anti-H_GIPR hIgG1 Antibody(AMG-133)	
ACVR2A:ACTRIIB:Active A	
ACVR2A KO HEK-293 Cell Line	ACVR2B KO HEK-293 Cell Line
Activin A Reporter Cell Line	BRE Reporter 293 Cell Line
H_ACVR2A Reporter Cell Line	H_ACVR2B Reporter Cell Line
H_ACVR2B Reporter DDX35TM Cell Line	H_ACVR2A HEK-293(ACVR2B KO) Cell Line
H_ACVR2B CHO-K1 Cell Line	H_ACVR2B HEK-293(ACVR2A KO) Cell Line
Anti-ACVR2A hIgG1 Antibody(LAE-102)	Anti-ACVR2B hIgG1 Antibody(Bimagrumab)
Anti-ACVR2B hIgG1 Antibody(Fab-17G05)	Anti-ACVR2B mIgG2a Antibody(Bimagrumab)
Anti-GDF8 hIgG4 Reference Antibody (Apitbio)	Anti-H_ACVR2B hIgG1 Reference Antibody(Bimbio)
Biotinylated Human ACVR2A Protein; His-Avi Tag	Biotinylated Human ACVR2B Protein; His-Avi Tag
Biotinylated Mouse ACVR2A Protein; His-Avi Tag	Biotinylated Mouse ACVR2B Protein; His-Avi Tag
Human Activin A Protein; His Tag	Human Activin A Protein; His Tag (CHO)
Human Activin B Protein; His Tag	Human ACVR2A Protein; hFc Tag
Human ACVR2A Protein; hFc Tag (Sotatercept)	Human ACVR2A Protein; His Tag
Human ACVR2B Protein; hFc Tag	Human ACVR2B Protein; His Tag
Human latent GDF-8 Protein; His Tag	Mouse ACVR2A Protein; His Tag
Mouse ACVR2B Protein; His Tag	
AMY:CALCR RAMP	
H_CALCR RAMP1(AMY1) Reporter CHO-K1 Cell Line	H_CALCR RAMP2(AMY2) Reporter CHO-K1 Cell Line
H_CALCR RAMP3(AMY3) Reporter CHO-K1 Cell Line	H_CALCR RAMP3(AMY3) β -Arrestin Reporter CHO-K1 Cell Line
H_CALCR Reporter CHO-K1 Cell Line	H_CALCR β -Arrestin Reporter CHO-K1 Cell Line
Rat_CALCR RAMP1(AMY1) Reporter COS-7 Cell Line	Rat_CALCR RAMP3(AMY3) Reporter COS-7 Cell Line
Rat_CALCR Reporter COS-7 Cell Line	
THRB	
H_THR β Reporter HEK-293 Cell Line	
MC4R	
H_MC4R Reporter HEK-293 Cell Line	
ASGR1	

H_ASGR1 CHO-K1 Cell Line	H_ASGR1 HEK-293 Cell Line
Anti-ASGR1 hIgG1 Antibody(4A2.001)	
Cynomolgus ASGR1 Protein; His Tag	Human ASGR1 Protein; hFc Tag
Human ASGR1 Protein; His Tag	Mouse ASGR1 Protein; His Tag

License Agreement:

By purchasing and using this cell line product, the user voluntarily agrees to accept and abide by the following policies:

- This cell line product is restricted to research use only and shall not be used for any commercial purposes.
- This product is strictly prohibited from being used in the diagnosis or treatment of human or animal diseases, and shall not be directly used in experiments involving humans.
- Users and their contractors engaged for their benefit may use this material and its derivatives only within the agreed research scope; modification of the material is not permitted, nor may it be distributed, sold, transferred, or otherwise provided to any other entity (including affiliates).
- If use beyond the above scope is required, prior written permission from Genomeditech (Shanghai) Co.,Ltd. must be obtained. For details, please contact Genomeditech (Shanghai) Co.,Ltd.