

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Product Sheet

H_FcRn CHO-K1 Cell Line

Catalog number: GM-C09637

Version 3.3.1.250526

H_FcRn CHO-K1 Cell Line is a clonal stable CHO-K1 cell line that at constitutively

expresses human FcRn and human B2M genes, constructed using lentiviral technology.

Quantity 5E6 Cells per vial,1 mL

Product Format 1 vial of frozen cells

Shipping Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

Target Human_FcRn & Human_B2M

Gene ID/Uniprot ID P55899 & P61769

Host Cell CHO-K1

Recovery Medium F12K+10% FBS+1% P.S

Growth medium F12K+10% FBS+1% P.S+100 μg/mL Hygromycin+4 μg/mL Puromycin

Note None

Freezing Medium 90% FBS+10% DMSO

Growth properties Adherent

Growth Conditions 37°C, 5% CO₂

Mycoplasma Testing The cell line has been screened to confirm the absence of Mycoplasma species.

Safety considerations Biosafety Level 2

Note It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.



Order: +86 021-68455258/50432826/50432825


Toll-free: +86 400 627 9288

Email: service@genomeditech.com

Materials

Reagent	Manufacturer/Catalogue No.
F12K	BOSTER/PYG0036
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Hygromycin	Genomeditech/GM-040403
Puromycin	Genomeditech/GM-040401
Anti-FcRn hIgG4 Reference Antibody(Rozabio)	Genomeditech/GM-87351MAB
Abflo® 488 Rabbit anti-Human β2 Microglobulin mAb	abclonal/A24223

Figures

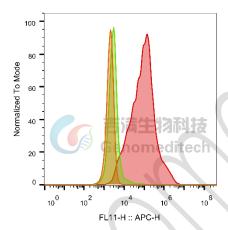

Sampleid	Geometric Mean : FLTT-H
CHO-K1 anti-FcRn+APC-2nd Ab	1969
CHO-K1 H_FcRn H_IgG+APC-2nd Ab	2914
CHO-K1 H_FcRn anti-FcRn+APC-2nd Ab	81228

Figure 1 | H_FcRn CHO-K1 Cell Line (Cat. GM-C09637) was determined by flow cytometry using Anti-FcRn hIgG4 Reference Antibody (Rozabio) (Cat. GM-87351MAB).



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com



10⁴ FL6-H :: FITC-H

SampleID	Geometric Mean : FL6-H
CHO-K1 anti-B2M+FITC-2nd Ab	1855
CHO-K1 H_FcRn Rabbit_lgG+FITC-2nd Ab	2665
CHO-K1 H_FcRn anti-B2M+FITC-2nd Ab	18996
	•

Figure 2 | H_FcRn CHO-K1 Cell Line (Cat. GM-C09637) was determined by flow cytometry using anti-Human β 2 Microglobulin mAb (abclonal/A24223).

10

Cell Recovery

Recovery Medium: F12K+10% FBS+1% P.S

10

102

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- e) Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Cell passage

Growth medium: F12K+10% FBS+1% P.S+100 μg/mL Hygromycin+4 μg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Remove and discard culture medium.
- b) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- c) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes at 37°C).
- d) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- e) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- f) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- g) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:4 - 1:5 is recommended

Medium Renewal: Every 2 to 3 days

Notes

a) After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

Sequence

FcRn(FCGRT) P55899

MGVPRPQPWALGLLIFLLPGSLGAESHLSLLYHLTAVSSPAPGTPAFWVSGWLGPQQYLSYNSLRGEAEPCG AWVWENQVSWYWEKETTDLRIKEKLFLEAFKALGGKGPYTLQGLLGCELGPDNTSVPTAKFALNGEEFMN FDLKQGTWGGDWPEALAISQRWQQQDKAANKELTFLLFSCPHRLREHLERGRGNLEWKEPPSMRLKARPSS PGFSVLTCSAFSFYPPELQLRFLRNGLAAGTGQGDFGPNSDGSFHASSSLTVKSGDEHHYCCIVQHAGLAQPL RVELESPAKSSVLVVGIVIGVLLLTAAAVGGALLWRRMRSGLPAPWISLRGDDTGVLLPTPGEAQDADLKDV NVIPATA

B2M P61769

 $MSRSVALAVLALLSLSGLEAIQRTPKIQVYSRHPAENGKSNFLNCYVSGFHPSDIEVDLLKNGERIEKVEHSDL\\ SFSKDWSFYLLYYTEFTPTEKDEYACRVNHVTLSQPKIVKWDRDM$

Related Products

FcγR	
Cynomolgus_FcRn MDCK Cell Line	H_FCGR1A(CD64) CHO-K1 Cell Line
H_FCGR1A(CD64) HEK-293 Cell Line	H_FCGR2A(CD32A) CHO-K1 Cell Line
H_FCGR2B(CD32B) CHO-K1 Cell Line	H_FCGR3A(CD16a) 158F CHO-K1 Cell Line



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

H_FCGR3A(CD16a) 158V CHO-K1 Cell Line	H_FCGR3B(CD16b) CHO-K1 Cell Line	
Mouse_FcRn MDCK Cell Line		
Anti-FcRn hIgG4 Reference Antibody(Rozabio)	Anti-H_FcRn IgG4 Antibody(Rozanolixizumab)	
Anti-Mouse CD1632 mIgG2b Antibody(2.4G2)		
ADCCP		
ADCC FcγRIIIa(158F) Jurkat Effector Cell Line	ADCC FcγRIIIa(158V) DDX35TM Jurkat Effector Cell Line	
ADCC FcγRIIIa(158V) Jurkat Effector Cell Line	ADCC M_FcγRIV Jurkat Effector Cell Line	
ADCP FcγRI Jurkat Effector Cell Line	ADCP FcγRIIa DDX35TM Jurkat Effector Cell Line	
ADCP FcγRIIa Jurkat Effector Cell Line	ADCP FcγRIIa R131 Jurkat Effector Cell Line	
ADCP FcγRIIb Jurkat Effector Cell Line		

Limited Use License Agreement

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.