

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

# **Product Sheet**


# Tango-H\_LPAR1 CHO-K1 Cell Line

Catalog number: GM-C14343

Version 3.3.1.251114

The LPAR1 gene (Lysophosphatidic Acid Receptor 1) encodes a G protein-coupled receptor (GPCR) that serves as a receptor for lysophosphatidic acid (LPA). It mediates several signaling pathways, such as Rho, Ras, PLC, PI3K, and MAPK, regulating cell proliferation, migration, survival, differentiation, and cytoskeletal remodeling. LPAR1 is widely expressed, particularly in the nervous, cardiovascular, and renal systems, and plays key roles in embryonic and neural development, fibrosis, and tumor metastasis, making it an important molecule in signal transduction and a potential therapeutic target.

Tango-H\_LPAR1 CHO-K1 Cell Line is a clonal stable CHO-K1 cell line constructed using lentiviral technology, constitutive expression of the LPAR1 gene constructed through Tango technology, along with signal-dependent expression of a luciferase reporter gene. When igand binds to LPAR1, it activates downstream signaling pathways, leading to the expression of luciferase. The luciferase activity measurement indicates the activation level of the signaling pathway and can thus be used to evaluate the in vitro effects of related drugs.

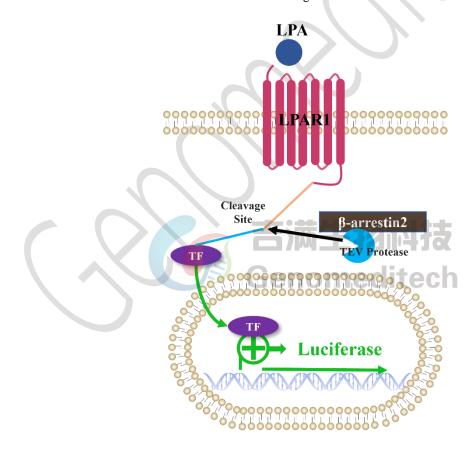



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

### **Specifications**

**Quantity** 5E6 Cells per vial,1 mL

**Product Format** 1 vial of frozen cells

**Shipping** Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

**Recovery Medium** F12K+10% FBS+1% P.S

Puromycin

Note None

Freezing Medium 90% FBS+10% DMSO

**Growth properties** Adherent

**Growth Conditions** 37°C, 5% CO<sub>2</sub>

Mycoplasma Testing The cell line has been screened to confirm the absence of Mycoplasma species.

**Safety considerations** Biosafety Level 2

Note It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.

#### **Materials**

Reagent	Manufacturer/Catalogue No.
F12K	BOSTER/PYG0036
Fetal Bovine Serum	ExCell/FSP500
Pen/Strep	Thermo/15140-122
Blasticidin	Genomeditech/GM-040404
Hygromycin	Genomeditech/GM-040403
Puromycin	Genomeditech/GM-040401
1-Oleoyl lysophosphatidic acid sodium	MCE/HY-107614
Anti-LPAR1 hIgG1 Antibody(h309-4-A-BP1)	Genomeditech/GM-87978AB
GMOne-Step 2.0 Luciferase Reporter Gene Assay Kit	Genomeditech/GM-040513

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

### **Figures**



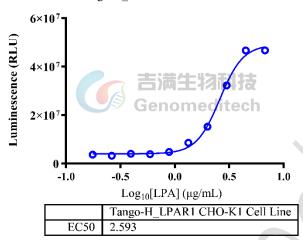


Figure 1 | Response to 1-Oleoyl lysophosphatidic acid sodium. The Tango-H\_LPAR1 CHO-K1 Cell Line (Cat. GM-C14343) at a concentration of 1E4 cells/well (96-well format) was stimulated with serial dilutions of 1-Oleoyl lysophosphatidic acid sodium(LPA) (MCE/HY-107614) in assay buffer (F12K+1% P.S) for 16 hours. The firefly luciferase activity was measured using the Luciferase Reporter Assay Kit (Genomeditech). The maximum induction fold was approximately [14.5]. Data are shown by drug mass concentration.


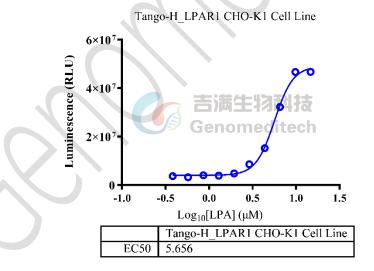


Figure 2 | Response to 1-Oleoyl lysophosphatidic acid sodium. The Tango-H\_LPAR1 CHO-K1 Cell Line (Cat. GM-C14343) at a concentration of 1E4 cells/well (96-well format) was stimulated with serial dilutions of 1-Oleoyl lysophosphatidic acid sodium(LPA) (MCE/HY-107614) in assay buffer (F12K + 1% FBS + 1% P.S) for 16 hours. The firefly luciferase activity was measured using the Luciferase Reporter Assay Kit (Genomeditech). The maximum induction fold was approximately [14.5]. Data are shown by drug molar concentration.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

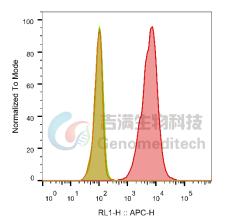

Г	SampleID	Geometric Mean : RL1-H
	CHO Anti-LPAR1+APC-2nd	93.6
	Tango-H-LPAR1 CHO H_lgG+APC-2nd	97.2
	Tango-H-LPAR1 CHO Anti-LPAR1+APC-2nd	5711

Figure 3 | Tango-H\_LPAR1 CHO-K1 Cell Line (Cat. GM-C14343) was determined by flow cytometry using Anti-LPAR1 hIgG1 Antibody(h309-4-A-BP1) (Cat. GM-87978AB).

# **Cell Recovery**

Recovery Medium: F12K+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- e) Incubate the culture at 37°C in a suitable incubator. A 5% CO<sub>2</sub> in air atmosphere is recommended if using the medium described on this product sheet.

## **Cell Freezing**

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

吉满生物科技 Genomeditech Genomeditech (Shanghai) Co.,Ltd.

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

#### Cell passage

Growth medium: F12K+10% FBS+1% P.S+4 µg/mL Blasticidin+100 µg/mL Hygromycin+4 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Remove and discard culture medium.
- b) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- c) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes at 37°C).
- d) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- e) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- f) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- g) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:4 - 1:5 is recommended

Medium Renewal: Every 2 to 3 days

#### **Notes**

a) After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

#### **Related Products**

LPAR1	
H_LPAR1 CHO-K1 Cell Line	

### **License Agreement:**

By purchasing and using this cell line product, the user voluntarily agrees to accept and abide by the following policies:

- This cell line product is restricted to research use only and shall not be used for any commercial purposes.
- This product is strictly prohibited from being used in the diagnosis or treatment of human or animal diseases, and shall not be directly
  used in experiments involving humans.
- Users and their contractors engaged for their benefit may use this material and its derivatives only within the agreed research scope; modification of the material is not permitted, nor may it be distributed, sold, transferred, or otherwise provided to any other entity (including affiliates).



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

If use beyond the above scope is required, prior written permission from Genomeditech (Shanghai) Co.,Ltd. must be obtained. For details, please contact Genomeditech (Shanghai) Co.,Ltd.