

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Product Sheet


H_ACVR2A Reporter Cell Line

Catalog number: GM-C28074

Version 3.3.1.250914

Activin type II receptors (ActRII) include two distinct receptors, ActRIIA (ACVR2A) and ActRIIB (ACVR2B), which are members of the transforming growth factor-beta (TGF- β) receptor family. They serve as type II receptors for other TGF- β family members, including Activin A, Activin B, BMP7, BMP9, BMP10, GDF1, GDF8 (myostatin), GDF11, and Nodal. ALK4 and ALK7 are the primary type I TGF- β receptor family members for Activin A and Activin B, respectively. Bimagrumab is a human monoclonal antibody that can block activin type II receptors (ActRII), thereby regulating skeletal muscle growth. Studies have shown that Bimagrumab has an affinity for both ActRIIB and ActRIIA, but its affinity for ActRIIB is at least 200 times stronger than for ActRIIA.

H_ACVR2A Reporter Cell Line is a clonal stable cell line that knockout ACVR2B and constitutive expression of human ACVR2A gene, along with signal-dependent expression of a luciferase reporter gene. When Activin A or GDF-8 binds to ACVR2A, it activates downstream signaling pathways, leading to the expression of luciferase. Blockade antibodies can inhibit this signal transmission. The luciferase activity measurement indicates the activation level of the signaling pathway and can thus be used to evaluate the in vitro effects of drugs related to this signaling pathway.

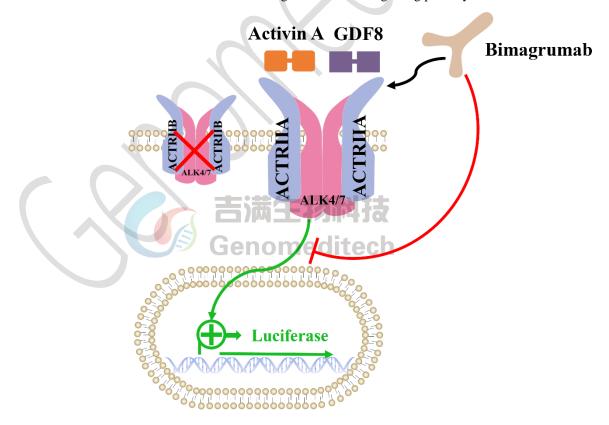



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

Specifications

Quantity 5E6 Cells per vial,1 mL

Product Format 1 vial of frozen cells

Shipping Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

Recovery Medium DMEM+10% FBS+1% P.S

Growth medium DMEM+10% FBS+1% P.S+400 μg/mL G418+125 μg/mL Hygromycin

These cells are constructed based on ACVR2B knockout parent cells, which contain Note

Blasticidin and Puromycin resistance genes.

Freezing Medium 90% FBS+10% DMSO

Growth properties Adherent

Growth Conditions 37°C, 5% CO₂

Mycoplasma Testing The cell line has been screened to confirm the absence of Mycoplasma species.

Safety considerations Biosafety Level 2

Note It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.

Materials

Reagent	Manufacturer/Catalogue No.
DMEM	Gibco/C11995500BT
Fetal Bovine Serum	ExCell/FSP500
Pen/Strep	Thermo/15140-122
G418	Genomeditech/GM-040402
Hygromycin	Genomeditech/GM-040403
Recombinant Human/Mouse/Rat GDF-8	Novoprotein/CJ43
Human Activin A Protein; His Tag	Genomeditech/GM-87616RP
Anti-ACVR2B hIgG1 Antibody(Fab-17G05)	Genomeditech/GM-82350AB
Anti-ACVR2B hIgG1 Antibody(Bimagrumab)	Genomeditech/GM-51148AB
GMOne-Step 2.0 Luciferase Reporter Gene Assay Kit	Genomeditech/GM-040513



Email: service@genomeditech.com

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Figures

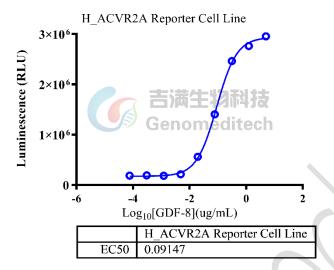


Figure 1 | Response to Recombinant Human/Mouse/Rat GDF-8. H_ACVR2A Reporter Cell Line (Cat. GM-C28074) at a concentration of 1.5E4 cells/well (96-well format) was stimulated with serial dilutions of Recombinant Human/Mouse/Rat GDF-8 (Novoprotein/CJ43) in assay buffer (DMEM + 1% FBS + 1% P.S) for 7 hours. The firefly luciferase activity was measured using the Luciferase Reporter Assay Kit (Genomeditech). The maximum induction fold was approximately [16.3]. Data are shown by drug mass concentration.

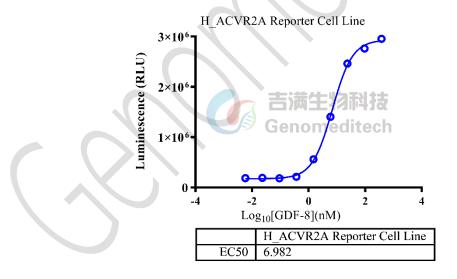


Figure 2 | Response to Recombinant Human/Mouse/Rat GDF-8. H_ACVR2A Reporter Cell Line (Cat. GM-C28074) at a concentration of 1.5E4 cells/well (96-well format) was stimulated with serial dilutions of Recombinant Human/Mouse/Rat GDF-8 (Novoprotein/CJ43) in assay buffer (DMEM + 1% FBS + 1% P.S) for 7 hours. The firefly luciferase activity was measured using the Luciferase Reporter Assay Kit (Genomeditech). The maximum induction fold was approximately [16.3]. Data are shown by drug molar concentration.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

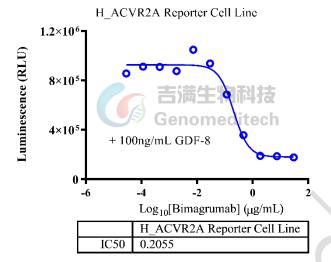


Figure 3 | Response to Anti-ACVR2B hIgG1 Antibody (Bimagrumab). Serial dilutions of antibodies were incubated with 1.5E4 cells/well of the H_ACVR2A Reporter Cell Line (Cat. GM-C28074) in a 96-well plate for 1 hour. Then, 100 ng/mL of Recombinant Human/Mouse/Rat GDF-8 (Novoprotein/CJ43) was added, and the incubation continued in assay buffer (DMEM + 1% FBS + 1% P.S) for 6 hours. Firefly luciferase activity is measured using the Luciferase Reporter Assay Kit (Genomeditech). The results indicated a maximum blocking fold of approximately [4.5]. Data are shown by drug mass concentration.

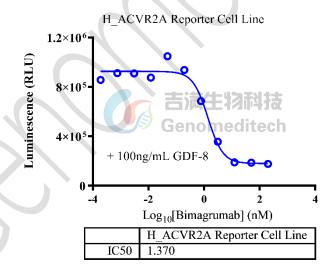


Figure 4 | Response to Anti-ACVR2B hIgG1 Antibody (Bimagrumab). Serial dilutions of antibodies were incubated with 1.5E4 cells/well of the H_ACVR2A Reporter Cell Line (Cat. GM-C28074) in a 96-well plate for 1 hour. Then, 100 ng/mL of Recombinant Human/Mouse/Rat GDF-8 (Novoprotein/CJ43) was added, and the incubation continued in assay buffer (DMEM + 1% FBS + 1% P.S) for 6 hours. Firefly luciferase activity is measured using the Luciferase Reporter Assay Kit (Genomeditech). The results indicated a maximum blocking fold of approximately [4.5]. Data are shown by drug molar concentration.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288



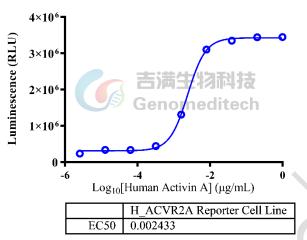


Figure 5 | Response to Human Activin A Protein. H_ACVR2A Reporter Cell Line (Cat. GM-C28074) at a concentration of 1.5E4 cells/well (96-well format) was stimulated with serial dilutions of Human Activin A Protein (Cat. GM-87616RP) in assay buffer (DMEM + 1% FBS + 1% P.S) for 6 hours. The firefly luciferase activity was measured using the Luciferase Reporter Assay Kit (Genomeditech). The maximum induction fold was approximately [12.7]. Data are shown by drug mass concentration.

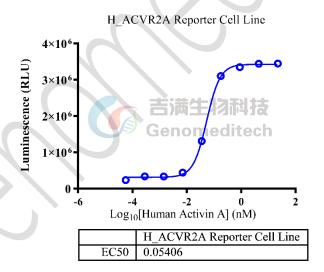


Figure 6 | Response to Human Activin A Protein. H_ACVR2A Reporter Cell Line (Cat. GM-C28074) at a concentration of 1.5E4 cells/well (96-well format) was stimulated with serial dilutions of Human Activin A Protein (Cat. GM-87616RP) in assay buffer (DMEM + 1% FBS + 1% P.S) for 6 hours. The firefly luciferase activity was measured using the Luciferase Reporter Assay Kit (Genomeditech). The maximum induction fold was approximately [12.7]. Data are shown by drug molar concentration.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

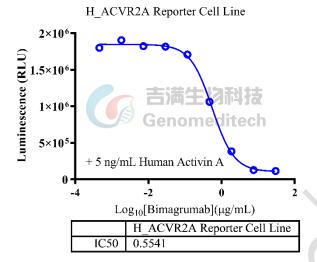


Figure 7 | Response to Anti-ACVR2B hIgG1 Antibody (Bimagrumab). Serial dilutions of antibodies were incubated with 1.5E4 cells/well of the H_ACVR2A Reporter Cell Line (Cat. GM-C28074) in a 96-well plate for 1 hour. Then, 5 ng/mL of Human Activin A Protein (Cat. GM-87616RP) was added, and the incubation continued in assay buffer (DMEM + 1% FBS + 1% P.S) for 6 hours. Firefly luciferase activity is measured using the Luciferase Reporter Assay Kit (Genomeditech). The results indicated a maximum blocking fold of approximately [16.2]. Data are shown by drug mass concentration.


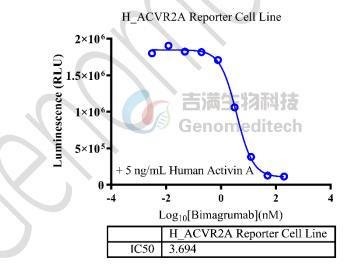


Figure 8 | Response to Anti-ACVR2B hIgG1 Antibody (Bimagrumab). Serial dilutions of antibodies were incubated with 1.5E4 cells/well of the H_ACVR2A Reporter Cell Line (Cat. GM-C28074) in a 96-well plate for 1 hour. Then, 5 ng/mL of Human Activin A Protein (Cat. GM-87616RP) was added, and the incubation continued in assay buffer (DMEM + 1% FBS + 1% P.S) for 6 hours. Firefly luciferase activity is measured using the Luciferase Reporter Assay Kit (Genomeditech). The results indicated a maximum blocking fold of approximately [16.2]. Data are shown by drug molar concentration.

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

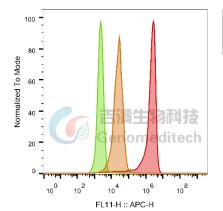

SampleID	Geometric Mean : FL11-H
Null anti-ACVR2B+APC+2nd Ab	23401
H_ACVR2A Reporter H_lgG+APC+2nd Ab	2274
H_ACVR2A Reporter anti-ACVR2B+APC+2nd Ab	1.48E6

Figure 9 | H_ACVR2A Reporter Cell Line (Cat. GM-C28074) was determined by flow cytometry using Anti-ACVR2B hIgG1 Antibody(Bimagrumab) (Cat. GM-51148AB).

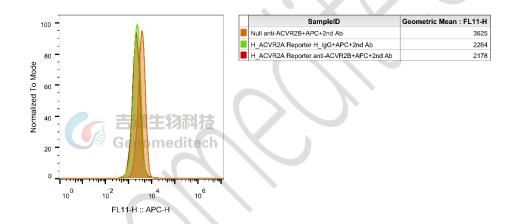


Figure 10 | H_ACVR2A Reporter Cell Line (Cat. GM-C28074) was determined by flow cytometry using Anti-ACVR2B hIgG1 Antibody(Fab-17G05) (Cat. GM-82350AB).

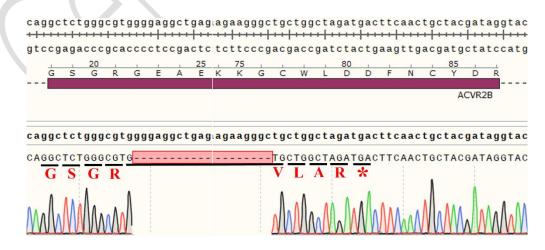


Figure 11 | The Sanger sequencing of the H_ACVR2A Reporter Cell Line showed successful knockout of ACVR2B.

吉満生物科技 Genomeditech

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

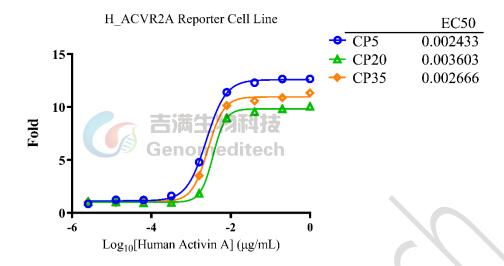


Figure 12 | The passage stability of response to Human Activin A Protein. The passage 5, 20 and 35 of H_ACVR2A Reporter Cell Line (Cat. GM-C28074) at a concentration of 1.5E4 cells/well (96-well format) were stimulated with serial dilutions of Human Activin A Protein (Cat. GM-87616RP) in assay buffer (DMEM + 1% FBS + 1% P.S) for 6 hours. The firefly luciferase activity was measured using the Luciferase Reporter Assay Kit (Genomeditech). Data are shown by drug mass concentration.

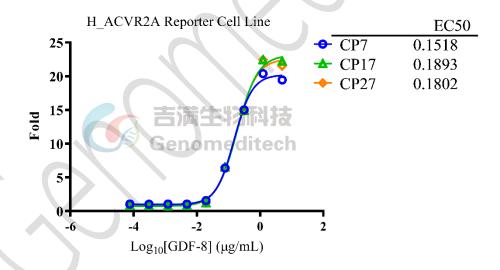


Figure 13 | The passage stability of response to Recombinant Human/Mouse/Rat GDF-8. The passage 7, 17 and 27 of H_ACVR2A Reporter Cell Line (Cat. GM-C28074) at a concentration of 1.5E4 cells/well (96-well format) were stimulated with serial dilutions of Recombinant Human/Mouse/Rat GDF-8 (Novoprotein/CJ43) in assay buffer (DMEM + 1% FBS + 1% P.S) for 6 hours. The firefly luciferase activity was measured using the Luciferase Reporter Assay Kit (Genomeditech). Data are shown by drug mass concentration.

吉满生物科技 Genomeditech Genomeditech (Shanghai) Co.,Ltd.

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

Cell Recovery

Recovery Medium: DMEM+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- e) Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: DMEM+10% FBS+1% P.S+400 µg/mL G418+125 µg/mL Hygromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Subculturing is necessary when the cell density reaches 80%. It is recommended to perform subculturing at a ratio of 1:3 to 1:4 every 2-3 days. Ensure that the density does not exceed 80%, as overcrowding can lead to reduced viability due to compression.
- b) Remove and discard culture medium.
- c) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- d) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 30 to 60 seconds at 37°C).
- e) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- f) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

g) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.

h) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:4 is recommended

Medium Renewal: Every 2 to 3 days

Notes

a) Upon initial thawing, a higher number of dead cells is observed, which is a normal phenomenon. Significant improvement is seen after adaptation. Once the cells reach a stable state, the number of dead cells decreases after subculturing and the cell growth rate becomes stable.

b) Ensure that the cell density does not exceed 80%, as overcrowding may lead to reduced viability due to compression.

Related Products

GCGR			
H_GCGR Reporter CHO-K1 Cell Line	H_GCGR Reporter HEK-293 Cell Line		
H_GCGR Reporter HEK-293 DDX35TM Cell Line	Cynomolgus_GCGR HEK-293 Cell Line		
H_GCGR CHO-K1 Cell Line	H_GCGR HEK-293 Cell Line		
Mouse_GCGR HEK-293 Cell Line			
Anti-H_GCGR hIgG2 Antibody(volagidemab)			
GLP1R			
H_GLP1R Reporter CHO-K1 Cell Line	H_GLP1R Reporter HEK-293 Cell Line		
H_GLP1R Reporter HEK-293 DDX35TM Cell Line	H_GLP1R β-Arrestin Reporter CHO-K1 Cell Line		
Cynomolgus_GLP1R HEK-293 Cell Line	H_GLP1R CHO-K1 Cell Line		
H_GLP1R GIPR CHO-K1 Cell Line	H_GLP1R HEK-293 Cell Line		
Mouse_GLP1R HEK-293 Cell Line			
Anti-GLP1R hIgG1 Antibody(mAb-36986)	Anti-H_GLP1R hIgG1 Antibody(glutazumab)		
FGFR1			
H_FGF21 Reporter HEK-293 Cell Line			
Human FGF-21 Protein; His Tag			
CALCA(CGRP): CALCRL RAMP			
H_CALCRL RAMP1 Reporter HEK-293 Cell Line	H_CALCRL RAMP1 Reporter HEK-293 DDX35TM Cell Line		
Cynomolgus_CALCRL RAMP1 HEK-293 Cell Line	H_CALCRL RAMP1 CHO-K1 Cell Line		
H_CALCRL RAMP1 HEK-293 Cell Line			
Anti-CALCRL RAMP1 hIgG2 Antibody(Erenumab)			
GIP:GIPR			
H_GIPR Reporter CHO-K1 Cell Line	H_GIPR Reporter HEK-293 Cell Line		
H_GIPR Reporter HEK-293 DDX35TM Cell Line	Cynomolgus_GIPR CHO-K1 Cell Line		
Cynomolgus_GIPR HEK-293 Cell Line	H_GIPR CHO-K1 Cell Line		
H_GIPR HEK-293 Cell Line	Mouse_GIPR CHO-K1 Cell Line		



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

Mouse_GIPR HEK-293 Cell Line			
Anti-H_GIPR hIgG1 Antibody(AMG-133)			
ACVR2A: ACTRIIB: Active A			
ACVR2A KO HEK-293 Cell Line	Activin A Reporter Cell Line		
BRE Reporter 293 Cell Line	H_ACVR2B Reporter Cell Line		
ACVR2B KO HEK-293 Cell Line	H_ACVR2A HEK-293(ACVR2B KO) Cell Line		
H_ACVR2B CHO-K1 Cell Line	H_ACVR2B HEK-293(ACVR2A KO) Cell Line		
	Anti-ACVR2B hIgG1 Antibody(Bimagrumab)		
Anti-ACVR2B hIgG1 Antibody(Fab-17G05)	Anti-ACVR2B mIgG2a Antibody(Bimagrumab)		
Anti-H_ACVR2B hIgG1 Reference Antibody(Bimbio)			
Biotinylated Human ACVR2A Protein; His-Avi Tag	Biotinylated Human ACVR2B Protein; His-Avi Tag		
Biotinylated Mouse ACVR2A Protein; His-Avi Tag	Biotinylated Mouse ACVR2B Protein; His-Avi Tag		
Human Activin A Protein; His Tag	Human Activin A Protein; His Tag (CHO)		
Human Activin B Protein; His Tag	Human ACVR2A Protein; hFc Tag		
Human ACVR2A Protein; hFc Tag (Sotatercept)	Human ACVR2A Protein; His Tag		
Human ACVR2B Protein; hFc Tag	Human ACVR2B Protein; His Tag		
Human latent GDF-8 Protein; His Tag	Mouse ACVR2A Protein; His Tag		
Mouse ACVR2B Protein; His Tag			
AMY: CALCR RAMP			
H_CALCR RAMP3(AMY3) Reporter CHO-K1 Cell Line	H_CALCR Reporter CHO-K1 Cell Line		

License Agreement:

By purchasing and using this cell line product, the user voluntarily agrees to accept and abide by the following policies:

- This cell line product is restricted to research use only and shall not be used for any commercial purposes.
- This product is strictly prohibited from being used in the diagnosis or treatment of human or animal diseases, and shall not be directly
 used in experiments involving humans.
- Users and their contractors engaged for their benefit may use this material and its derivatives only within the agreed research scope; modification of the material is not permitted, nor may it be distributed, sold, transferred, or otherwise provided to any other entity (including affiliates).
- If use beyond the above scope is required, prior written permission from Genomeditech (Shanghai) Co.,Ltd. must be obtained. For details, please contact Genomeditech (Shanghai) Co.,Ltd.