

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Product Sheet

H_CD40L CHO-K1 Cell Line

Catalog number: GM-C35532

Version 3.3.1.250605

H_CD40L CHO-K1 Cell Line is a clonal stable CHO-K1 cell line that constitutively **Description**

expresses the human CD40 gene, constructed using lentiviral technology.

Quantity 5E6 Cells per vial,1 mL

Product Format 1 vial of frozen cells

Shipping Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

Target Human_CD40L

Gene ID/Uniprot ID P29965

Host Cell CHO-K1

Recovery Medium F12K+10% FBS+1% P.S

Growth medium F12K+10% FBS+1% P.S+4 μg/mL Puromycin

Note None

Freezing Medium 90% FBS+10% DMSO

Growth properties Adherent

Growth Conditions 37°C, 5% CO₂

Mycoplasma Testing The cell line has been screened to confirm the absence of Mycoplasma species.

Safety considerations Biosafety Level 2

Note It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Materials

| Reagent | Manufacturer/Catalogue No. |
|--|-----------------------------|
| F12K | BOSTER/PYG0036 |
| Fetal Bovine Serum | Cegrogen biotech/A0500-3010 |
| Pen/Strep | Thermo/15140-122 |
| Puromycin | Genomeditech/GM-040401 |
| Anti-CD40L hIgG1 Reference Antibody (Frebio) | Genomeditech/GM-87982MAB |
| Anti-H_CD40L hIgG1 Antibody(dapirolizumab) | Genomeditech/GM-83977AB |

Figures

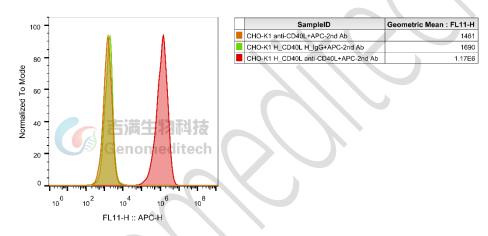


Figure 1 | H_CD40L CHO-K1 Cell Line (Cat. GM-C35532) was determined by flow cytometry using Anti-CD40L hIgG1 Reference Antibody (Frebio) (Cat. GM-87982MAB).

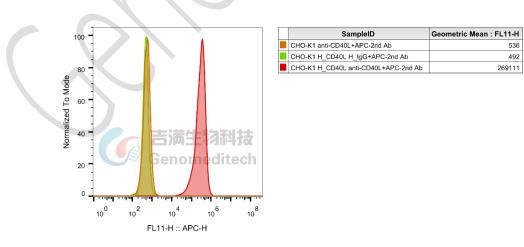



Figure 2 | H_CD40L CHO-K1 Cell Line (Cat. GM-C35532) was determined by flow cytometry using Anti-H_CD40L hIgG1 Antibody(dapirolizum) (Cat. GM-83977AB).



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

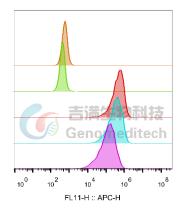

| Geometric Mean : FL11-H |
|-------------------------|
| 572 |
| 440 |
| 4.17E5 |
| 304762 |
| 111264 |
| |

Figure 3 | The passage stability of the H_CD40L CH0-K1 Cell Line (Cat. GM-C35532) was determined by flow cytometry using Anti-H_CD40L hIgG1 Antibody(dapirolizum) (Cat. GM-83977AB).

Cell Recovery

Recovery Medium: F12K+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- e) Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

Cell passage

Growth medium: F12K+10% FBS+1% P.S+4 μ g/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Remove and discard culture medium.
- b) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- c) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes at 37°C).
- d) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- e) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- f) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- g) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:4 - 1:5 is recommended

Medium Renewal: Every 2 to 3 days

Notes

a) After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

Sequence

CD40L(CD40LG) P29965

MIETYNQTSPRSAATGLPISMKIFMYLLTVFLITQMIGSALFAVYLHRRLDKIEDERNLHEDFVFMKTIQRCNT GERSLSLLNCEEIKSQFEGFVKDIMLNKEETKKENSFEMQKGDQNPQIAAHVISEASSKTTSVLQWAEKGYYT MSNNLVTLENGKQLTVKRQGLYYIYAQVTFCSNREASSQAPFIASLCLKSPGRFERILLRAANTHSSAKPCGQ QSIHLGGVFELQPGASVFVNVTDPSQVSHGTGFTSFGLLKL

Related Products

| CD40: CD40L | |
|--|---|
| H_CD40(TNFRSF5) Reporter 293 Cell Line | H_CD40L HEK-293 Cell Line |
| H_CD40(TNFRSF5) Reporter Jurkat Cell Line | Cynomolgus_CD40 CHO-K1 Cell Line |
| Cynomolgus_CD40L CHO-K1 Cell Line | H_CD40(TNFRSF5) CHO-K1 Cell Line |
| H_CD40(TNFRSF5) HEK-293 Cell Line | |
| Anti-H_CD40 hIgG1 Antibody(APX005M) | Anti-H_CD40 hIgG1 Antibody(ravagalimab) |
| Anti-H_CD40L hIgG1 Antibody(dapirolizumab) | Anti-H_CD40L hIgG1 Antibody(frexalimab) |
| Biotinylated Human CD40 Protein; His-Avi Tag | Cynomolgus CD40 Protein; His Tag |
| Human CD40 Protein; His Tag | Human CD40L Protein; His Tag |



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

| IFN-α | | |
|--|--------------------------------------|--|
| IFNα Reporter HEK-293 Cell Line | IFNα Reporter MDCK Cell Line | |
| IFNα Reporter THP1 Cell Line | | |
| BCMA:BAFFR:TACI | | |
| H_BAFFR Jurkat Blockade Reporter Cell Line | H_BAFFR Reporter Cell Line | |
| H_BCMA Reporter Cell Line | H_TACI Reporter Cell Line | |
| Cynomolgus_BCMA CHO-K1 Cell Line | H_BCMA CHO-K1 Cell Line | |
| H_BCMA HEK-293 Cell Line | | |
| Anti-BAFF hIgG1 Antibody(belimumab) | Anti-BAFFR hIgG1 Antibody(ianalumab) | |
| Anti-BCMA hIgG1 Antibody(Belantamab) | Anti-BCMA hIgG1 Antibody(SEA-BCMA) | |
| Anti-BCMA hIgG4 Antibody(BCMB69) | | |
| Biotinylated Human BAFF Protein; His-Avi Tag | Cynomolgus BAFF Protein; His Tag | |
| Human BAFF Protein; His Tag | Mouse BAFF Protein; His Tag | |
| BDCA2(CLEC4C) | | |
| H_BDCA2 Jurkat Cell Line | H_BDCA2 Reporter Jurkat Cell Line | |
| Cynomolgus_BDCA2 CHO-K1 Cell Line | Cynomolgus_BDCA2 Jurkat Cell Line | |
| H_BDCA2 CHO-K1 Cell Line | H_BDCA2 HEK-293 Cell Line | |
| Anti-H_BDCA2 hIgG1 Antibody(Litifilimab) | | |
| Cynomolgus BDCA2 Protein; His Tag | Human BDCA2 Protein; His Tag | |

Limited Use License Agreement

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.