

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Product Sheet

Cynomolgus_IL-23R HEK-293 Cell Line

Catalog number: GM-C39197

Version 3.3.1.250221

Cynomolgus_IL-23R HEK-293 Cell Line is a clonal stable HEK-293 cell line that

Description constitutively expresses the cynomolgus IL-23R gene, constructed using lentiviral

technology.

Quantity 5E6 Cells per vial,1 mL

Product Format 1 vial of frozen cells

Shipping Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

Target Cynomolgus_IL-23R

Gene ID/Uniprot ID G7NWY5

Host Cell HEK-293

Recovery Medium DMEM+10% FBS+1% P.S

Growth medium DMEM+10% FBS+1% P.S+0.75 μg/mL Puromycin

Note None

Freezing Medium 90% FBS+10% DMSO

Growth properties Adherent

Growth Conditions 37°C, 5% CO₂

Mycoplasma Testing The cell line has been screened to confirm the absence of Mycoplasma species.

Safety considerations Biosafety Level 2

Note It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Materials

Reagent	Manufacturer/Catalogue No.
DMEM	Gibco/C11995500BT
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Puromycin	Genomeditech/GM-040401
Anti-IL-23R hIgG1 Antibody(5D4)	Genomeditech/GM-87786AB

Figures

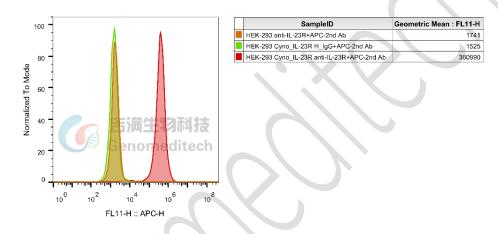


Figure 1 | Cynomolgus_IL-23R HEK-293 Cell Line (Cat. GM-C39197) was determined by flow cytometry using Anti-IL-23R hIgG1 Antibody(5D4) (Cat. GM-87786AB).

Cell Recovery

Recovery Medium: DMEM+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.

吉满生物科技 Genomeditech Genomeditech (Shanghai) Co.,Ltd.

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

e) Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: DMEM+10% FBS+1% P.S+0.75 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Subculturing is necessary when the cell density reaches 80%. It is recommended to perform subculturing at a ratio of 1:3 to 1:4 every 2-3 days. Ensure that the density does not exceed 80%, as overcrowding can lead to reduced viability due to compression.
- b) Remove and discard culture medium.
- c) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- d) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 30 to 60 seconds at 37°C).
- e) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- f) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- g) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- h) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:4 is recommended

Medium Renewal: Every 2 to 3 days

Notes

- a) Upon initial thawing, a higher number of dead cells is observed, which is a normal phenomenon. Significant improvement is seen after adaptation. Once the cells reach a stable state, the number of dead cells decreases after subculturing and the cell growth rate becomes stable.
- b) Ensure that the cell density does not exceed 80%, as overcrowding may lead to reduced viability due to compression.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Sequence

IL-23R(IL23R) G7NWY5

MNQVTIQWDAVIALYILFSWCHGGITNINCSGHIWVEPATIFKMGMNISIYCQAAIKNCQPRQLHFYKNGIKE RFEITRINKTTARLWYKNFLEPHASMYCTAECPNHFQETLICGKDISSGYPPDVPDEVTCVIYEYSGNMTCTW NAGKLTYVDTKYVVHVKSLETEEEQQYLTSSYINISTDSLQGGKKYLVWVQAENALGMEESKQLQIHLDDIV IPSASIISRAETINATVPKTIIYWDSQTTIEKVSCEMRYKATTNQTWNVKEFDTNFTYVQQSEFYLEPNIKYIFQ VRCQETGKRYWQAWSSPFFHKTTETVPQVTSKSSQHNTWNSGLTVASISTGHLTSDNREDIGLLLGMIVFAV MLSILSLIGIFNRSLRTGIKRRILLLIPKWLYEDIPNMKNSNVVKMLQENNELMNNNSSEQVLYVDPTITEIKEIF IPEHKPTDYKKENTGPLETRDYLQNSLFDNTTVVYIPDLNTGYKPQISNFLPGRSRLSNNDEITSLTLKPPVDSL DSGNNPRSQKHPNFAFSVSSVNSPSSTLFLGELSLILNQGECSSPDIQNSVEEETTMLLENDSPSETIPEQTLLPD EFVSCLGIMNEELPSINSYFPQSILESHFKRISLLEK

Related Products

IL-23		
H_IL-23 Reporter 293 Cell Line	H_IL-23R HEK-293 Cell Line	
Human IL-23A&IL-12B Heterodimer Protein; His Tag	Human IL-23R Protein; hFc Tag	
TNF:TNFR2:TNFR1		
H_TNFR2 Null Reporter Cell Line	H_TNFR2 Reporter Jurkat Cell Line	
H_TNFR2 Reporter V2 Cell Line	Cynomolgus_TNFRSF1B(TNFR2) CHO-K1 Cell Line	
H_TNFRSF1B(TNFR2) CHO-K1 Cell Line	H_TNFRSF1B(TNFR2) HEK-293 Cell Line	
Membrane Bound H_TNFα CHO-K1 Cell Line	Membrane Bound H_TNFα(cleavage-resistant) CHO-K1 Cell Line	
Anti-H_TNFR2 hIgG1 Antibody(1H10)	Anti-H_TNFRSF1B(TNFR2) hIgG1 Antibody(UC2.3.8)	
Anti-TNFR1 hIgG1 Antibody(Atrosab)	Anti-TNF- a hIgG1 Antibody (CT-P17)	
TL1A:DR3(TNFRSF25)		
H_TNFRSF25(DR3) Reporter Jurkat Cell Line	H_TNFSF15(TL1A) Reporter Cell Line	
Mouse_TNFRSF25(DR3) Reporter Jurkat Cell Line	Cynomolgus_TNFSF15(TL1A) HEK-293 Cell Line	
H_TNFRSF25(DR3) CHO-K1 Cell Line	H_TNFRSF25(DR3) HEK-293 Cell Line	
H_TNFSF15(TL1A) CHO-K1 Cell Line	H_TNFSF15(TL1A) HEK-293 Cell Line	
Mouse_TNFSF15(TL1A) HEK-293 Cell Line		
Anti-H_TNFSF15(TL1A) hIgG1 Antibody(PF-06480605)	Anti-H_TNFSF15(TL1A) hIgG1 Antibody(Tulisokibart、PRA-023)	
Anti-H_TNFSF15(TL1A) hIgG4 Antibody	Anti-TL1A hIgG1 Reference Antibody (Duvbio)	
Anti-TL1A hIgG1 Reference Antibody (Tulbio)	Anti-TNFRSF25(DR3) hIgG1 Antibody(PTX-35)	
Biotinylated Cynomolgus TL1A Protein; His-Avi Tag	Cynomolgus TL1A Protein; His Tag	
Human TL1A Protein; His Tag		



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

Limited Use License Agreement

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.