

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

# **Product Sheet**

## H\_IL33R HEK-293 Cell Line

Catalog number: GM-C39950

Version 3.3.1.250624

H\_IL33R HEK-293 Cell Line is a clonal stable HEK-293 line that constitutively expresses

the human IL33R gene, constructed using lentiviral technology.

**Quantity** 5E6 Cells per vial,1 mL

**Product Format** 1 vial of frozen cells

**Shipping** Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

Target H\_IL33R

Gene ID/Uniprot ID Q01638-1

Host Cell HEK-293

**Recovery Medium** DMEM+10% FBS+1% P.S

Growth medium DMEM+10% FBS+1% P.S+0.75 μg/mL Puromycin

Note None

Freezing Medium 90% FBS+10% DMSO

**Growth properties** Adherent

Growth Conditions 37°C, 5% CO<sub>2</sub>

**Mycoplasma Testing** The cell line has been screened to confirm the absence of Mycoplasma species.

**Safety considerations** Biosafety Level 2

Note It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

#### **Materials**

Reagent	Manufacturer/Catalogue No.
DMEM	Gibco/C11995500BT
Fetal Bovine Serum	ExCell/FSP500
Pen/Strep	Thermo/15140-122
Puromycin	Genomeditech/GM-040401
IL-33R (ST2) Monoclonal Antibody (hIL33Rcap),	APC, Thermofisher/17-9338-42
eBioscience <sup>TM</sup>	

## **Figures**

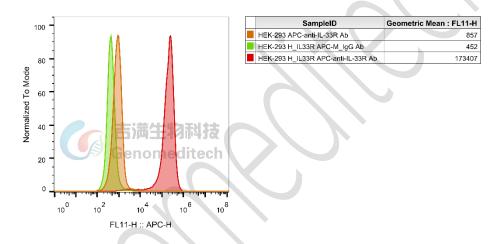


Figure 1 | H\_IL33R HEK-293 Cell Line (Cat. GM-C39950) was determined by flow cytometry using IL-33R (ST2) Monoclonal Antibody (hIL33Rcap), APC, eBioscience™ (Thermofisher/17-9338-42).

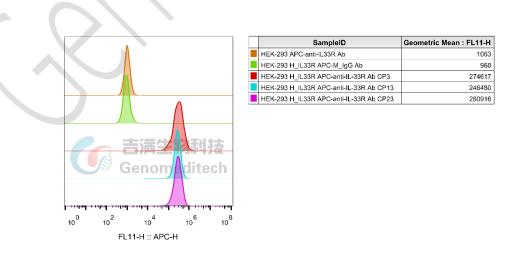


Figure 2 | The passage stability of the H\_IL33R HEK-293 Cell Line (Cat. GM-C39950) was determined by flow cytometry using IL-33R (ST2) Monoclonal Antibody (hIL33Rcap), APC, eBioscience™ (Thermofisher/17-9338-42).

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

**Cell Recovery** 

Recovery Medium: DMEM+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival,

continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C.

Storage at -70°C will result in loss of viability.

Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring

and cap out of the water. Thawing should be rapid (approximately 2 - 3 minutes).

b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying

with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.

c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately

176 x g for 5 minutes. Discard supernatant.

d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.

e) Incubate the culture at 37°C in a suitable incubator. A 5% CO<sub>2</sub> in air atmosphere is recommended if using the medium

described on this product sheet.

**Cell Freezing** 

Freezing Medium: 90% FBS+10% DMSO

Centrifuge at 176 x g for 3 minutes to collect cells.

Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL. b)

Aliquot 1 mL into each vial. c)

d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid

nitrogen as soon as possible.

Cell passage

Growth medium: DMEM+10% FBS+1% P.S+0.75 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth

medium.

Subculturing is necessary when the cell density reaches 80%. It is recommended to perform subculturing at a ratio of

1:3 to 1:4 every 2-3 days. Ensure that the density does not exceed 80%, as overcrowding can lead to reduced viability

due to compression.

Remove and discard culture medium. b)

Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor. c)

d) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell

layer is dispersed (usually within 30 to 60 seconds at 37°C).

Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. e)

Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

- f) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- g) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- h) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:4 is recommended

Medium Renewal: Every 2 to 3 days

#### **Notes**

- a) Upon initial thawing, a higher number of dead cells is observed, which is a normal phenomenon. Significant improvement is seen after adaptation. Once the cells reach a stable state, the number of dead cells decreases after subculturing and the cell growth rate becomes stable.
- b) Ensure that the cell density does not exceed 80%, as overcrowding may lead to reduced viability due to compression.

## **Sequence**

IL33R(IL1RL1) Q01638-1

MGFWILAILTILMYSTAAKFSKQSWGLENEALIVRCPRQGKPSYTVDWYYSQTNKSIPTQERNRVFASGQLL KFLPAAVADSGIYTCIVRSPTFNRTGYANVTIYKKQSDCNVPDYLMYSTVSGSEKNSKIYCPTIDLYNWTAPL EWFKNCQALQGSRYRAHKSFLVIDNVMTEDAGDYTCKFIHNENGANYSVTATRSFTVKDEQGFSLFPVIGAP AQNEIKEVEIGKNANLTCSACFGKGTQFLAAVLWQLNGTKITDFGEPRIQQEEGQNQSFSNGLACLDMVLRI ADVKEEDLLLQYDCLALNLHGLRRHTVRLSRKNPIDHHSIYCIIAVCSVFLMLINVLVIILKMFWIEATLLWRD IAKPYKTRNDGKLYDAYVVYPRNYKSSTDGASRVEHFVHQILPDVLENKCGYTLCIYGRDMLPGEDVVTAV ETNIRKSRRHIFILTPQITHNKEFAYEQEVALHCALIQNDAKVILIEMEALSELDMLQAEALQDSLQHLMKVQ GTIKWREDHIANKRSLNSKFWKHVRYQMPVPSKIPRKASSLTPLAAQKQ\*

### **Related Products**

IL-33	
H_IL33 Reporter 293 Cell Line	H_IL33R CHO-K1 Cell Line
Anti-IL33 hIgG4 Antibody(Itepekimab)	Anti-IL33 hIgG4 Reference Antibody (Itepbio)

## **License Agreement:**

By purchasing and using this cell line product, the user voluntarily agrees to accept and abide by the following policies:

- This cell line product is restricted to research use only and shall not be used for any commercial purposes.
- This product is strictly prohibited from being used in the diagnosis or treatment of human or animal diseases, and shall not
  be directly used in experiments involving humans.
- Users are not permitted to modify the cell line in any way, nor to share, distribute, sell, sublicense, or otherwise transfer
  the licensed materials or their progeny to other laboratories, departments, research institutes, hospitals, universities,
  biotechnology companies, or any other third parties, except for research activities outsourced on behalf of the licensee.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

If the product is intended to be transferred to a third party, used for commercial development, preclinical or clinical drug functional validation, commercial production testing, or any other applications beyond the scope of this statement, prior written permission must be obtained from Genomeditech (Shanghai) Co.,Ltd. For details, please contact Genomeditech (Shanghai) Co.,Ltd.