

Genomeditech (Shanghai) Co.,Ltd.

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Product Sheet

Cynomolgus_CD8A CD8B CHO-K1 Cell Line

Catalog number: GM-C42780

Version 3.3.1.251114

Cynomolgus_CD8A CD8B CHO-K1 Cell Line is a clonal stable CHO-K1 cell line that

Description constitutively expresses the cynomolgus CD8A and cynomolgus CD8B gene, constructed

using lentiviral technology.

Quantity 5E6 Cells per vial,1 mL

Product Format 1 vial of frozen cells

Shipping Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

Target Cynomolgus_CD8A & Cynomolgus_CD8B

Gene ID/Uniprot ID XP_065382911.1 & XP_065382905.1

Host Cell CHO-K1

Recovery Medium F12K+10% FBS+1% P.S

Growth medium F12K+10% FBS+1% P.S+4 μg/mL Puromycin

Note None

Freezing Medium 90% FBS+10% DMSO

Growth properties Adherent

Growth Conditions 37°C, 5% CO₂

Mycoplasma Testing The cell line has been screened to confirm the absence of Mycoplasma species.

Safety considerations Biosafety Level 2

Note It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.



Genomeditech (Shanghai) Co.,Ltd.

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Materials

Reagent	Manufacturer/Catalogue No.
F12K	BOSTER/PYG0036
Fetal Bovine Serum	ExCell/FSP500
Pen/Strep	Thermo/15140-122
Puromycin	Genomeditech/GM-040401
Anti-CD8A hIgG4 Antibody(REGN5054)	Genomeditech/GM-88299AB
Anti-CD8B hIgG1 Antibody(xhCD8v12)	Genomeditech/GM-88302AB

Figures


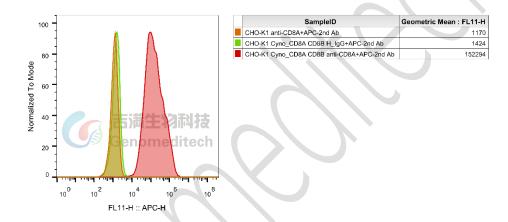


Figure 1 | Cynomolgus_CD8A CD8B CH0-K1 Cell Line (Cat. GM-C42780) was determined by flow cytometry using Anti-CD8A hIgG4 Antibody(REGN5054) (Cat. GM-88299AB).

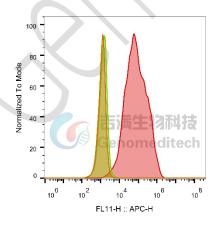

SampleID	Geometric Mean : FL11-H
CHO-K1 anti-CD8B+APC-2nd Ab	1168
CHO-K1 Cyno_CD8A CD8B H_IgG+APC-2nd Ab	1363
CHO-K1 Cyno_CD8A CD8B anti-CD8B+APC-2nd Ab	76563

Figure 2 | Cynomolgus_CD8A CD8B CHO-K1 Cell Line (Cat. GM-C42780) was determined by flow cytometry using Anti-CD8B hIgG1 Antibody(xhCD8v12) (Cat. GM-88302AB).

吉满生物科技 Genomeditech Genomeditech (Shanghai) Co.,Ltd.

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

Cell Recovery

Recovery Medium: F12K+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- e) Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: F12K+10% FBS+1% P.S+4 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Remove and discard culture medium.
- b) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- c) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes at 37°C).
- d) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- e) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- f) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- g) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:4 - 1:5 is recommended



Genomeditech (Shanghai) Co.,Ltd.

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

Medium Renewal: Every 2 to 3 days

Notes

a) After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

Sequence

CD8A XP 065382911.1

MAPPVTALLLPLVLLLHAARPSQFRVSPLGRTWNLGETVELKCQVLLSNPTSGCSWLFQPRGTAARPTFLLYL SQNKPKAAEGLDTQRFSGKRLGDTFVLTLRDFRQENEGYYFCSALSNSIMYFSRFVPVFLPAKPTTTPAPRPPT PAPTTASQPLSLRPEACRPAAGGSVNTRGLDFACDIYIWAPLAGACGVLLLSLVITLYCNHRNRKRVCKCPRT LFCCQGLE

CD8B XP_065382905.1

MRPRLWLLLAAQLAVVRGSSVLQQTPAYIQVQTNKMAMLSCEAKISLGNMRIYWLRQRQAPSSDSHHEFLA LWDSTKGTVHGEEVAQEKIAVFRDASRFILNLTSVKPEDSGIYFCMIIGSPELTFGKGTQLSVVDFLPTTAQPT KKSTPKKRGCRLPRPATQKGPLCSPITLGLLVAGVLVLLVSLGVAIHLYCELSPLGSSGVPVCCLKCSWDTYF SSGAEQLAAWWECASENLWPAEESPASFHQTASRGRYFRNPCPPMPARKLQQYCNLKEAA

Related Products

CD8		
H_CD8A CD8B HEK-293 Cell Line	H_CD8A HEK-293 Cell Line	
Anti-mouse CD8B(Lyt 3.2) RIgG1 Antibody(53-5.8)	Anti-mouse CD8α mIgG2a Antibody(53-6.7)	

License Agreement:

By purchasing and using this cell line product, the user voluntarily agrees to accept and abide by the following policies:

- This cell line product is restricted to research use only and shall not be used for any commercial purposes.
- This product is strictly prohibited from being used in the diagnosis or treatment of human or animal diseases, and shall not be directly
 used in experiments involving humans.
- Users and their contractors engaged for their benefit may use this material and its derivatives only within the agreed research scope;
 modification of the material is not permitted, nor may it be distributed, sold, transferred, or otherwise provided to any other entity (including affiliates).
- If use beyond the above scope is required, prior written permission from Genomeditech (Shanghai) Co.,Ltd. must be obtained. For details, please contact Genomeditech (Shanghai) Co.,Ltd.