

# Product Sheet

## H\_Insulin Receptor A(IRA) CHO-K1 Cell Line

Catalog number: GM-C43946

Version 3.3.1.260122

**Description** H\_Insulin Receptor A(IRA) CHO-K1 Cell Line is a clonal stable CHO-K1 cell line that constitutively expresses the H\_Insulin Receptor A(Isoform A) gene, constructed using lentiviral technology.

**Quantity** 5E6 Cells per vial, 1 mL

**Product Format** 1 vial of frozen cells

**Shipping** Shipped on dry ice

**Storage Conditions** Liquid nitrogen immediately upon receipt

**Target** Human Insulin Receptor A(IRA)

**Gene ID/Uniprot ID** P06213-2

**Host Cell** CHO-K1

**Recovery Medium** F12K+10% FBS+1% P.S

**Growth medium** F12K+10% FBS+1% P.S+4 µg/mL Puromycin

**Note** None

**Freezing Medium** 90% FBS+10% DMSO

**Growth properties** Adherent

**Growth Conditions** 37°C, 5% CO<sub>2</sub>

**Mycoplasma Testing** The cell line has been screened to confirm the absence of Mycoplasma species.

**Safety considerations** Biosafety Level 2

**Note** It is recommended to expand the cell culture and store a minimum of 10 vials at an early passage for potential future use.

## Materials

Reagent	Manufacturer/Catalogue No.
F12K	BOSTER/PYG0036
Fetal Bovine Serum	ExCell/FSP500
Pen/Strep	Thermo/15140-122
Puromycin	Genomeditech/GM-040401
Anti-H_INSR hIgG1 Antibody(Valanafusp)	Genomeditech/GM-27376AB

## Figures

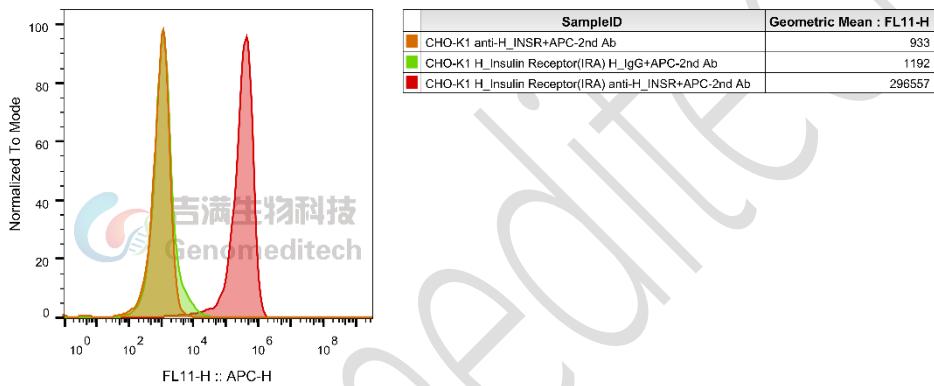


Figure 1 | H\_Insulin Receptor A(IRA) CHO-K1 Cell Line (Cat. GM-C43946) was determined by flow cytometry using Anti-H\_INSR hIgG1 Antibody(Valanafusp) (Cat. [GM-27376AB](#)).

## Cell Recovery

Recovery Medium: F12K+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 - 3 minutes).
- Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.

e) Incubate the culture at 37°C in a suitable incubator. A 5% CO<sub>2</sub> in air atmosphere is recommended if using the medium described on this product sheet.

## Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

## Cell passage

Growth medium: F12K+10% FBS+1% P.S+4 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Remove and discard culture medium.
- b) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- c) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes at 37°C).
- d) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- e) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- f) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- g) Incubate cultures at 37°C.

**Subcultivation Ratio: A subcultivation ratio of 1:4 - 1:5 is recommended**

**Medium Renewal: Every 2 to 3 days**

## Notes

- a) After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

## Sequence

H\_Insulin P06213-2

MATGGRRGAAAAPLLAVAALLGAAGHLYPGEVCPGMDIRNNLTRLHELENCSVIEGHLQILMFKTRPE  
DFRDLSFPKLIMITDYLLLFRVYGLESLKDLFPNLTVIRGSRLFFNYALVIFEMVHLKELGLYNLMNITRGSVRI  
EKNNELCYLATIDWSRILDSVEDNYIVLNKDDNEECGDCPGTAKGKTNCPATVINGQFVERCWTHSHCQKV  
CPTICKSHGCTAEGLCCHSECLGNCSQPDDPTKCVACRFYLDGRCVETCPPYYHFQDWRCVNFSCQDLH

HKCKNSRRQGCHQYVIHNNKCIPECPGTYMNSSNLLCTPCLGCPKVCHLLEGEKTIDSVTSQAQELRGCTVI  
 NGSLIINIRGGNNLAAELEANLGLIEEISGYLKIRRSYALVSLSSFRKLRLIRGETLEIGNYSFYALDNQNLRLQ  
 WDWSKHNLITQGKLFFHYNPKLCLSEIHKMEEVSGTGRQERNDIALKTNGDQASCENELLKFSYIRTSFD  
 KILLRWEPYWPPDFRDLLGFMLFYKEAPYQNVTEFDGQDACGSNSWTVVDIDPPLRSNDPKSQNHPGWL  
 GLKPWTQYAIFVKTIVTFSERRTYGAKSDIIYVQTDATNPSVPLDPISVNSSSQIILKWKP  
 PPSDPNGNITHYL  
 VFWERQAEDSELFELDYCLKGLKLPRTWSPPFESEDSQKH  
 NQSEYEDSAGECCSCP  
 KTDSQLKELEESSFRK  
 TFEDYLHNVVFVPRPSRKRRSLGDVGNVTAVPTVA  
 AFPNTSSTS  
 VPTSP  
 EEEHRP  
 FEKVV  
 NKESL  
 VISGL  
 RHFT  
 GYRIELQACNQDTPEERCSVA  
 AYVSARTMPEAKADDIVGPV  
 THEIFENN  
 VVHLM  
 WQEP  
 KEP  
 NGL  
 LIVL  
 YEVSY  
 RRYGDEELHLCVSRKHF  
 ALERGCRLRGLSPGN  
 SVRIRATSLAGNG  
 SWTEPTYF  
 YVTDYLD  
 VPSNIA  
 KIIIGPL  
 IFVFLFSV  
 VIGSIY  
 LFLRK  
 QPDG  
 PLG  
 PLY  
 YASSN  
 PEY  
 LSA  
 SDV  
 FPC  
 SVY  
 VP  
 DEWEV  
 SRE  
 KIT  
 LL  
 RELG  
 QGS  
 FGMV  
 YEGNARD  
 IIK  
 GEA  
 ETR  
 VAV  
 KTV  
 NES  
 ASLR  
 ERIE  
 FL  
 NEAS  
 VMKG  
 FTCH  
 H  
 V  
 RLL  
 G  
 V  
 VSK  
 G  
 Q  
 P  
 TL  
 V  
 M  
 E  
 L  
 M  
 A  
 H  
 GDL  
 KSY  
 L  
 R  
 S  
 L  
 R  
 P  
 E  
 A  
 E  
 N  
 N  
 P  
 G  
 R  
 P  
 P  
 T  
 L  
 Q  
 E  
 M  
 I  
 Q  
 M  
 A  
 A  
 E  
 I  
 A  
 D  
 G  
 M  
 A  
 Y  
 L  
 N  
 A  
 K  
 K  
 F  
 V  
 H  
 R  
 D  
 L  
 A  
 A  
 R  
 N  
 C  
 M  
 V  
 A  
 H  
 D  
 F  
 T  
 V  
 K  
 I  
 G  
 D  
 F  
 M  
 T  
 R  
 D  
 I  
 Y  
 E  
 T  
 D  
 Y  
 Y  
 R  
 K  
 G  
 G  
 K  
 L  
 P  
 V  
 R  
 W  
 M  
 A  
 P  
 E  
 S  
 L  
 K  
 D  
 G  
 V  
 F  
 T  
 T  
 S  
 S  
 D  
 M  
 W  
 S  
 F  
 G  
 V  
 V  
 L  
 W  
 E  
 I  
 T  
 S  
 L  
 A  
 E  
 Q  
 P  
 Y  
 Q  
 G  
 L  
 S  
 N  
 E  
 Q  
 V  
 L  
 K  
 F  
 V  
 M  
 D  
 G  
 G  
 Y  
 L  
 D  
 Q  
 P  
 D  
 N  
 C  
 P  
 E  
 R  
 V  
 T  
 D  
 L  
 M  
 R  
 M  
 C  
 W  
 Q  
 F  
 N  
 P  
 K  
 M  
 R  
 P  
 T  
 F  
 L  
 E  
 I  
 V  
 N  
 L  
 L  
 H  
 P  
 S  
 F  
 P  
 E  
 V  
 S  
 F  
 F  
 H  
 S  
 E  
 E  
 N  
 K  
 A  
 P  
 E  
 S  
 E  
 E  
 L  
 E  
 M  
 E  
 F  
 E  
 D  
 M  
 E  
 N  
 V  
 P  
 L  
 D  
 R  
 S  
 S  
 H  
 C  
 Q  
 R  
 E  
 E  
 A  
 G  
 G  
 R  
 D  
 G  
 G  
 S  
 S  
 L  
 G  
 F  
 K  
 R  
 S  
 Y  
 E  
 E  
 H  
 I  
 P  
 Y  
 T  
 H  
 M  
 N  
 G  
 G  
 K  
 K  
 N  
 G  
 R  
 I  
 L  
 T  
 L  
 P  
 R  
 S  
 N  
 P  
 S

## Related Products

INSR:insulin	
<a href="#">Insulin Receptor A(IRA) Reporter Cell Line</a>	<a href="#">Insulin Receptor B(IRB) Reporter Cell Line</a>
<a href="#">H_Insulin Receptor B(IRB) CHO-K1 Cell Line</a>	<a href="#">H_Insulin Receptor B(IRB) HEK-293 Cell Line</a>
<a href="#">Anti-H_INSR hIgG1 Antibody(Valanafusp)</a>	<a href="#">Anti-INSR hIgG1 Antibody(XPA.15.247.2.018)</a>

## License Agreement:

**By purchasing and using this cell line product, the user voluntarily agrees to accept and abide by the following policies:**

- This cell line product is restricted to research use only and shall not be used for any commercial purposes.
- This product is strictly prohibited from being used in the diagnosis or treatment of human or animal diseases, and shall not be directly used in experiments involving humans.
- Users and their contractors engaged for their benefit may use this material and its derivatives only within the agreed research scope; modification of the material is not permitted, nor may it be distributed, sold, transferred, or otherwise provided to any other entity (including affiliates).
- If use beyond the above scope is required, prior written permission from Genomeditech (Shanghai) Co.,Ltd. must be obtained. For details, please contact Genomeditech (Shanghai) Co.,Ltd.